

## Backing Up GASP Data and converting your data from version to version.

Backing up your data. - We recommend doing this prior to upgrading your software and periodically as it's possible for data to become corrupt or for hard drives to fail from time to time.

First of all, you need to know where your data is stored.

Usually data is stored in the below folders however this may differ slightly depending on where you initially installed your database.

GASPLab 3 data C:\ProgramData\InterVations\GASP Lab

GASPLab 4,5,6 data is stored in C:\ProgramData\InterVations\GASPLab4

\*\*Note Program data is a hidden file or folder and you will need to activate show hidden files or folders on your PC.

Windows 7 video <https://youtu.be/007vjyNqv8>

Windows 10 video <https://youtu.be/rBteQswtPI0>

However its better to double check the registry.

Go to start>run type regedit (windows 7) - Start type regedit (windows 10)

Open the registry editor and navigate to HKEY\_LOCAL\_MACHINE\SOFTWARE\WOW6432Node\InterVations\GaspLab4\Paths and note where the server, temp and video are stored

Server	REG_SZ	C:\ProgramData\InterVations\GASPLab4
Temp	REG_SZ	C:\ProgramData\InterVations\GASPLab4\Temp
Video	REG_SZ	C:\ProgramData\InterVations\GASPLab4\Video


Once you have located the folder copy the entire folder to another are od you PC or better to an external hard drive.

## Converting Data from GASPLab 3 to GASPLab 4

Install the latest version of GASPLab 3 [www.gaspsystems.com/download/\\_gasplab.exe](http://www.gaspsystems.com/download/_gasplab.exe)

Download and run [www.gaspsystems.com/download/convertgl3.exe](http://www.gaspsystems.com/download/convertgl3.exe)

Select the bottom option convert to mp4 and deinterlace and select Go this operation will take a long time to convert the data.

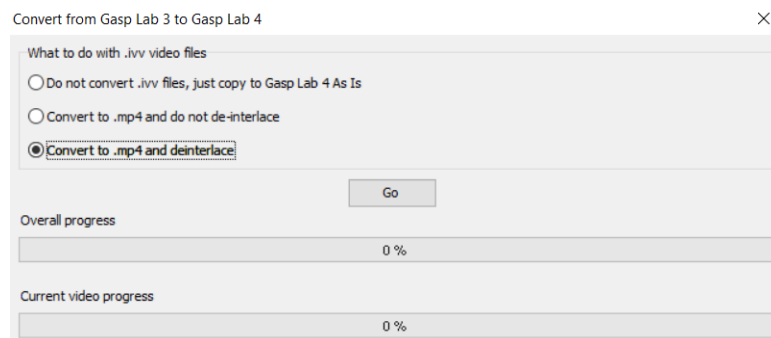

**Converting Data from GASPLab 4 to GASPLab 5**

**By installing GASPLab 5 the data will convert automatically. Please use the latest version**

[www.gaspsystems.com/download/\\_gasplab6.exe](http://www.gaspsystems.com/download/_gasplab6.exe)

Updating from GASPLab 5 to GASPLab 6 Please use the latest version

[www.gaspsystems.com/download/601/\\_gasplab6.exe](http://www.gaspsystems.com/download/601/_gasplab6.exe)